




# ALLERGARD<sup>ES</sup>



The AllerGard<sup>®</sup> NU-620 & NU-621 Animal Transfer Stations provide ISO Class 4 air protection and allergen protection through the use of vertical HEPA-filtered airflow on the work surface and an air barrier at the cabinet's opening(s). This dynamic airflow minimizes personnel exposure to animal hair/dander and cage bedding material.

## Standard Features:

- Single or Dual Access Options
- Two Washable or Disposable Single Mesh Pre-Filters
- 14" [356 mm] Hinged Access Opening, Easily Fits Rat Cages
- Motorized Auto-Lift Base Stand with Heavy-Duty Casters
- 150% HEPA Filter Loading Capacity
- AeroMax™ Control Center with Airflow Monitor & Alarm
- LED lighting with Red Light Option
- ISO Class 4 Air Protection

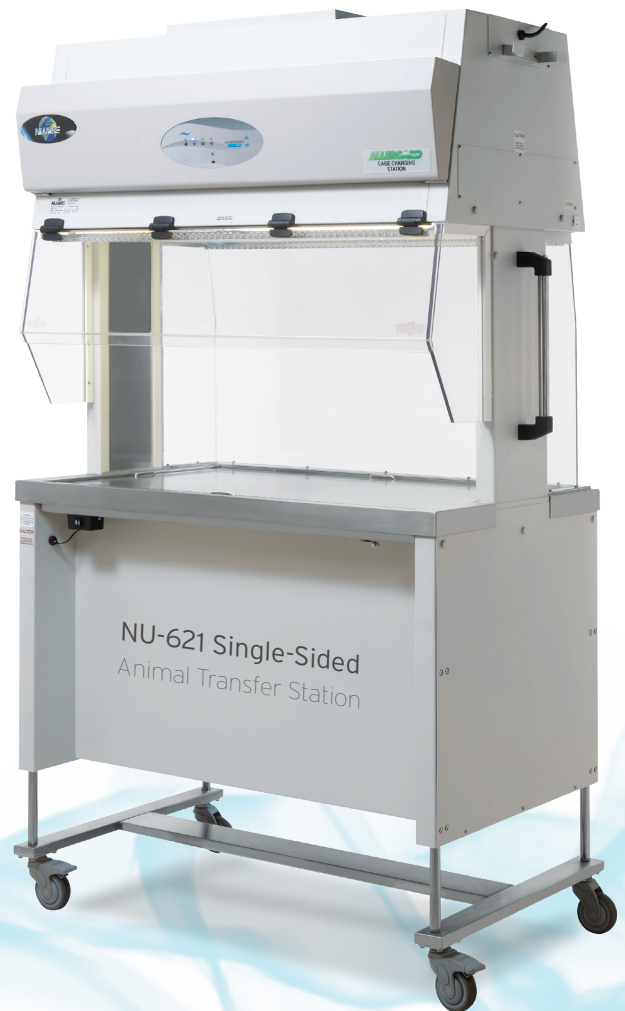

# NU-620 / NU-621 Specifications

Model #	NU-620	NU-621
Nominal Widths	3, 4, 5 ft [0.9, 1.2, 1.5 m]	
Work Zone Dimensions (W x D x H)	36.6 - 54.6 x 28 x 25.5 in [1067 - 1524 x 762 x 648 mm]	
Exterior Dimensions (W x D x H)	48 - 66 x 30.6 x 80.5 - 92.5 in [1220 - 1677 x 777 x 2045 - 2350 mm]	
Net Weight	430 - 530 lbs [195 - 240 kg]	
Certifications	UL listed, CE	
Access Opening / Type	Dual-Sided 14 in [356 mm] (Hinged)	Single-Sided 14 in [356 mm] (Hinged)
Electrical	115 VAC / 60 Hz   230 VAC / 50-60 Hz	
Filter Load Capacity	150%	
Airflow Dynamics	25 fpm [0.13 m/s] In-Flow, 60 fpm [0.30 m/s] Down-Flow	75 fpm [0.38 m/s] In-Flow, 60 fpm [0.30 m/s] Down-Flow
	ISO Class 4 Air Protection	
Supply Filter	3 in [76 mm] HEPA 99.99% Efficient @ 0.3 microns	
Pre-Filter	Two (2) Washable or Disposable Mesh Pre-Filter	
Exhaust Filter	3 in [76 mm] HEPA 99.99% Efficient @ 0.3 microns	
Control System	Aeromax™ w/ PressureFlow™ Pressure Sensor	
Lighting	LED	
	Red (custom)	
Construction Options	Removable Stainless Steel Prop-Up Work Tray	
	Side-Mounted Pull Bars	
	Stainless Steel Spill Trough	
	AttenuMont™ Vibration Control Support System	
	Electrical Outlet (Side-Wall Mounted)	
	Magnifying Glass Window	
	Mounted Hand Sanitizer Dispenser	
Integrated Base Stand	Stainless Steel Folding Shelf	
	Animal Feed Hoppers	
Integrated Base Stand	Motorized Adjustable Hydraulic Base Stand	

Standard Feature

Optional Feature

## Popular Features and Accessories




Removable Prop-Up Work Tray & Washable Mesh Pre-Filter



Stainless Steel Folding Shelf



Cage Lid Holder

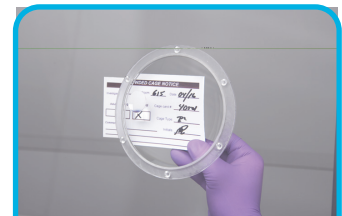

Window Sash Magnifying Glass



To Learn More or Request a Quote:

Visit: [www.nuaire.com](http://www.nuaire.com)

or Call: 1-763-553-1270

Toll Free: 1-800-328-3352

